


Allograft DBM Putty & DBM Crunch Putty

- 100% Allograft
 - Putty is cortical bone
 - Crunch Mix is cortical and cancellous bone
- Ready to use out of the package
- Resists irrigation
- All lots tested for BMP and osteoinductivity
- Demineralisation process validated for viral inactivation



Osteoinductive

Demineralised Bone Matrix (DBM) contains a variety of biologically active bone morphogenic proteins (BMPs), which are among the osteoinductive components of bone. More than 20 BMPs have been identified in human allograft bone tissue. A validated in vitro osteoinductive assay using ELISA BMP 2 concentration measurement method is used to verify that BMP is present in our final DBM product.

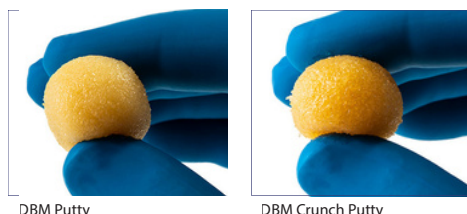

DBM Putty

DBM Crunch Putty

Product	Volume	Item Code
DBM Putty	0.5 cc	2014-40
	1.0 cc	2015-40
	2.5 cc	2016-40
	5.0 cc	2017-40
	10.0 cc	2018-40
DBM Crunch Putty	1.0 cc	4992-40
	2.5 cc	4993-40
	5.0 cc	4994-40
	10.0 cc	4995-40

Safety Assurance

Maxxeus DBM Putty and DBM Crunch Mix is packaged in a double, sterile pouch in a plunger. Each graft has a two year expiration date and is stored at ambient temperature. Grafts are prepared using a patented bioburden reduction process, which when combined with terminal irradiation provides a Sterility Assurance Level of 10⁻⁶. All donors are recovered and processed in the United States of America.

Speak to your local Business Development Manager for further information or contact us using the details below:

T: 01443 719 555

E: customerservice@hiuk.co.uk

www.hospitalinnovations.com

Hospital Innovations Limited

Concept House

Talbot Green Business Park

Pontyclun

CF72 9FG

